

mPIXdaq: Data acquisition for *miniPIX (EDU)* pixel detector

Vers. 0.9.3, September 2025

The miniPIX EDU is a camera for radiation based on the Timepix pixel read-out chip with 256x256 radiation-sensitive pixels of 55x55 μm^2 area and 300 μm depth each. The chip is covered by a very thin foil to permit α and β radiation to reach the pixels. The device is enclosed in an aluminum housing with a USB 2.0 interface. The sensor chip is covered by a thin foil and is very fragile; this area should be covered with a protective material if not measuring α radiation.

Other than single semi-conductor chips or simple Geiger counters, this device provides two-dimensional images of particle traces in the sensitive detector material. The high spatial resolution compared to the typical range of particles in silicon is useful to distinguish the different types of radiation and measure their deposited energies. α -particles are completely absorbed and deposit all of their energy in the sensitive area, allowing usage of the device as an energy spectrometer.

The vendor provides a ready-to-use program for different computer platforms as well as a software-development kit for own applications.

The code provided here is a minimalist example to read out single frames, i.e. a full set of 256x256 measurements accumulated over a given, fixed time interval. Each frame is displayed as an image with a logarithmic color scale representing the deposited energy. The analysis of the recorded signals, i.e. clustering of pixels, energy determination and visualization, is achieved with standard open-source tools for data analysis. It is therefore well-suited to give high-school or university students detailed insights and to enable them to carry out their own studies.

Getting ready for data taking

This code has been tested on *Ubuntu*, *openSuse*, *Fedora*, on Windows 64bit with *Python3.7.9* and on *Raspberry Pi* for the 32- and 64-bit versions of OS12. Other Linux distribution should not pose any problem.

On MS Windows, the libraries provided by the vendor only support *Python* vers. 3.7.9; such a rather historic version can be set up using e.g. the *miniconda* framework.

The code also supports devices other than the miniPIX EDU if the configuration files are available and copied to the *factory/* directory in the *pypixet Python* interface.

To get started, follow the steps below:

- Get the code from gitlab @ KIT

```
git clone https://gitlab.kit.edu/Guenter.Quast/mPIXdaq.
```

This repository includes the *Python* code and a minimalistic set of libraries provided by ADVACAM. `cd` to the *mPIXdaq* directory you just downloaded.

- Set up the USB interface of your computer to recognize the miniPIX EDU:

```
sudo install_driver_rules.sh
```

 (to be done only once), then connect the *miniPIX EDU* device to your computer.

The package may also be installed in your virtual python environment:

- `python -m pip install .`

Now everything is set up to enjoy your miniPIX EDU. Just run the *Python* program from any working directory by typing

```
run_mPIXdaq.py.
```

If you plan to record data, note that the path to the output file is relative to the current working directory.

Note that the *pypixet* initialization is set up to write log-files and configuration data to the directory */tmp/mPIX/*.

It is also worth mentioning that on some systems the current directory, ".", needs to be contained in the *LD_LIBRARY_PATH* so that the ADVACAM *Python* interface *pypixet* finds all its *C* libraries. This is also done in the *Python* script *run_mPIXdaq.py* by temporarily modifying the environment variable *LD_LIBRARY_PATH* if necessary and then restarting to execute the *Python* code in the new environment.

Starting the *Python* code by a different mechanism may not work without adjusting the environment variable *LD_LIBRARY_PATH*. In the *bash* shell, this is achieved by `export LD_LIBRARY_PATH='.'` on the command line. Note, however, that such a permanent change opens up a security gap on your computer!

Running the example script

Available options of the *Python* example to steer data taking and data archival to disk are shown by typing

`run_mPIXdaq.py --help`, resulting in the following output:

```
usage: run_mPIXdaq.py [-h] [-v VERBOSITY] [-o OVERLAY] [-a ACQ_TIME] [-c ACQ_COUNT] [-f FILE] [-w WRITEFILE]
                    [--circularity_cut CIRCULARITY_CUT] [-r READFILE]
```

read, analyze and display data from miniPIX device

options:

```
-h, --help                show this help message and exit
-v VERBOSITY, --verbosity VERBOSITY
                           verbosity level (1)
-o OVERLAY, --overlay OVERLAY
                           number of frames to overlay in graph (10)
-a ACQ_TIME, --acq_time ACQ_TIME
                           acquisition time/frame (0.1)
-c ACQ_COUNT, --acq_count ACQ_COUNT
                           number of frames to add (5)
-f FILE, --file FILE      file to store frame data
-w WRITEFILE, --writefile WRITEFILE
                           csv file to write cluster data
-t TIME, --time TIME      run time in seconds
--circularity_cut CIRCULARITY_CUT
                           circularity cut
-r READFILE, --readfile READFILE
                           file to read frame data
```

The default values are adjusted to situations with low rates, where frames from the *miniPIX* with an integration time of `acq_time = 0.5 s` are read. For the graphics display, `overlay = 10` recent frames are overlaid, leading to a total integration time of 5 s. These images represent a two-dimensional pixel map with a color code indicating the energy measured in each pixel.

Data analysis consists of clustering of pixels in each overlay-frame and determination of cluster parameters, like the number of pixels, energy and circularity. The threshold on circularity is controlled by the parameter `circularity_cut` ranging from 0. for perfectly linear to 1. for perfectly circular clusters. Technically, the covariance matrix of the clusters is calculated, and the circularity is defined as the ratio of the smaller and the larger of the two eigenvalues of the covariance matrix. This simple procedure already provides a good separation of α and β particles and of isolated pixels not assigned to clusters. The latter ones have a high probability of being produced in interactions of photons, while electrons from β radiation or from photon interactions produce long traces, and α particles produce large, circular clusters due to their very high ionization loss in the detector material.

A further, very sensitive variable is the variance of the energy distribution in the clusters. For α particles, the distribution peaks at the centre and steeply falls off towards the boundary, leading to a very small variance. A small ratio of the variances of the energy distribution and of the area covered by pixels is therefore a very prominent signature of α particles.

Properties of clusters are optionally written to a file in `.csv` format for later off-line analysis. A *Jupyter* notebook, *analyze_mPIXclusters.ipynb*, is provided which illustrates an example analysis.

To test the software without access to a miniPIX device or without a radioactive source, a file with recorded data is provided. Use the option `--readfile data/BlackForestStone.npy.gz` to start a demonstration. Note that the analysis of the recorded pixel frames is done in real time and may take some time on slow computers.

Implementation Details

The default data acquisition is based on the function *doSimpleIntegralAcquisition()* from the *ADVACAM Python* API. A fixed number of frames (`acq_counts`) with an adjustable accumulation time (`acq_time`) are read from the miniPIX device and added up.

The chosen readout mode is *PX_TPXMODE_TOT*, where “ToT” means “time over threshold”. This quantity shows good proportionality to the deposited energy at high signal values, but exhibits a strong non-linear behavior for small signals near the detection threshold of the miniPIX. Calibration constants are stored on the miniPIX device for each pixel, which are used to provide deposited energies per pixel in units of keV.

The relevant libraries for device control are provided in directories *advacam_<arch>* for x86_64 Linux, arm32 and arm64 and for Macintosh arm64 and MS Windows architectures. The contents of a typical directory is:

```
__init__.py    # package initialization
pypixet.so     # the Pixet Python interface
minipix.so     # C library for pypixet
pxcore.so      # C library for pypixet
pixet.ini      # initialization file, in same directory as pypixet
factory/       # initialization constants
```

Note that the copyright of these libraries belongs to ADVACAM. The libraries may be downloaded from their web page, ADVACAM DWONLOADS. They are provided here as *Python* packages for some platforms for convenience.

Data Analysis

The analysis shown in this example is intentionally very simple and based on standard libraries and functions. Pixels clustering is performed by finding connected regions in the pixel image with *scipy.ndimage.label()*. The shape of the clusters is determined from the ratio of the smaller and the larger one of the two eigenvalues of the covariance matrix calculated from the *x* and *y* coordinates of the pixels in a cluster. For circular clusters, as typically produced by α radiation, this ratio is close to one, while it is almost zero for the longer traces from β radiation.

The figure below shows the graphical display of the program with a pixel image and the typical distributions of the pixel and cluster energies and the number of pixels per cluster. The source used was a weakly radioactive stone from the Black Forest containing a small amount of Uranium and its decay products. The pixel map shown in the figure was sampled over a time of five seconds. The histogram in the lower-right corner shows how well the cluster types discriminate different types of radiation: α rays in the green band with relatively low numbers of pixels per cluster, electrons (β) as long tracks with large numbers of pixels per cluster and rather low energies. Single pixels not associated to clusters originate from γ rays. Some of the electron tracks with typically low energies also stem from photon interactions in the detector material (via the Compton process).

The analysis shown here is suitable for low-rate scenarios, e.g. the analysis of natural radiation as emitted by minerals like Pitchblend (=Uraninit), Columbit, Thorianit and others. Radon accumulated from the air in basement rooms on the surface of an electrostatically charged ballon also work fine. Therefore, the frame collection is chosen to be on the order of seconds, so that analysis results can be displayed in real-time on a sufficiently fast computer including the Raspberry Pi 5.

For applications at higher rates, the analysis may have to be done off-line by reading data from recorded files, or multiple cores must be used for the analysis task.

Sensor Details

The miniPIX (EDU) is based on the *Timepix* hybrid silicon pixel device, consisting of a semiconductor detector chip segmented into 256 x 256 square pixels with a pitch of 55 mm that is bump-bonded to the readout chip. Each element of the pixel matrix is connected to its own preamplifier, discriminator and digital counter integrated on the readout chip.

The built-in *Medipix2* variant of the chip is operated in the so-called “frame mode”, i.e. all pixels are read out at the same time, providing one frame consisting of the deposited energies per pixel collected during the acquisition time. If operated in time-over-threshold (ToT) mode, returned pixel readings represent the time the signal is over a given threshold in counts of the chip clock (appr. 10 MHz). *ToT* is linearly related to the energy deposition for large deposits exceeding 50 keV. The functional dependence on the deposited energy *E*, including threshold effects, is approximated by the following function

$$ToT = a E + b - c / (E - t)$$

Approximate values of the calibrations constants are *a* = 1.6, *b*=23, *c*=23 and *t*=4.3. Each pixel has its individual calibration stored on the chip, which is optionally applied to obtain pixel readings in units of keV. The calibration is reliable up to pixel energies of one MeV. Higher pixel energies may result when frames with short acquisition time

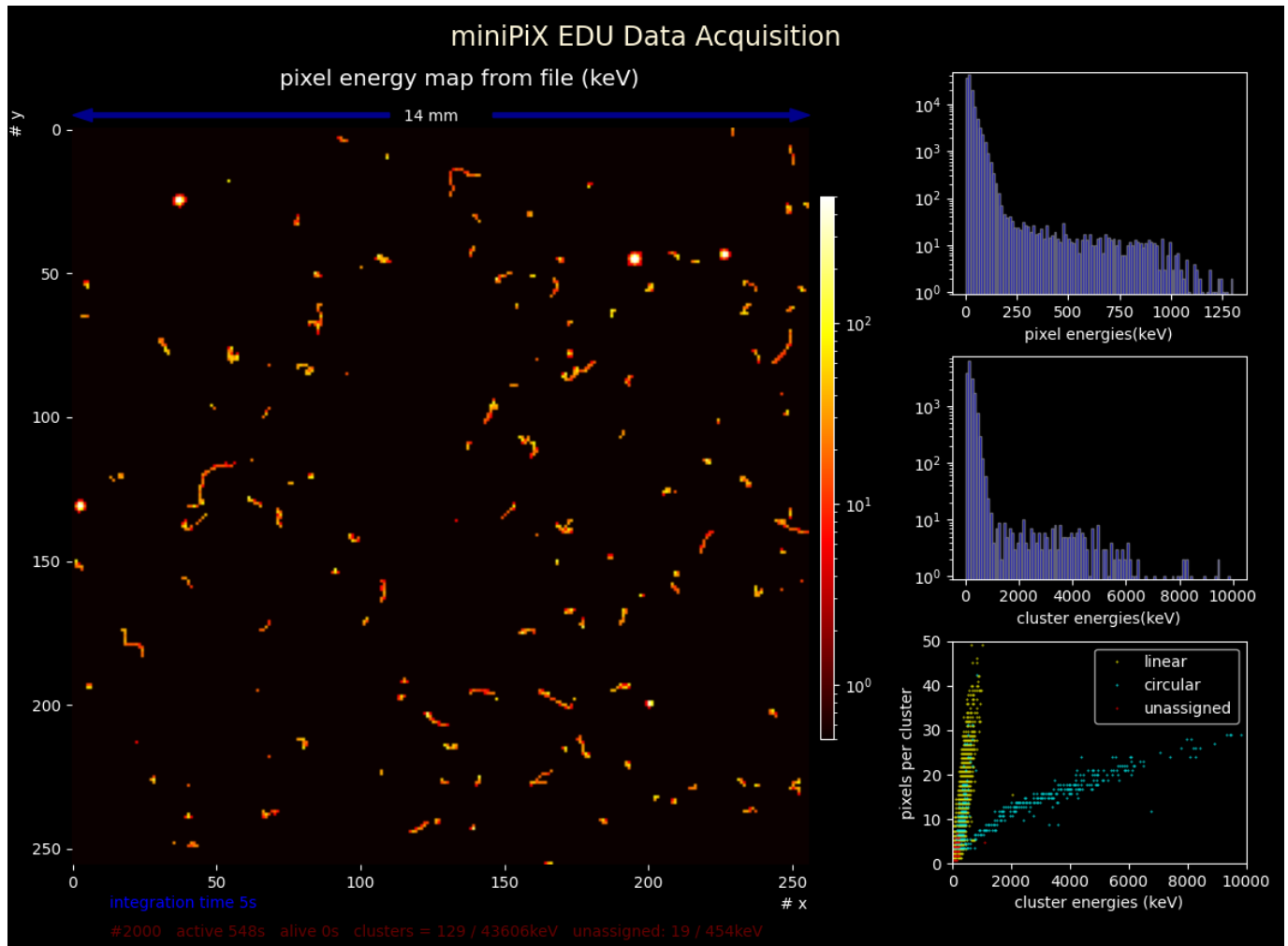


Abbildung 1: The graphical display of miniPiXdaq

are summed up. For details, see the article by J. Jakubek, *Precise energy calibration of pixel detector working in time-over-threshold mode*, NIM A 633 (2011), 5262-5265*.

Package Structure

This package consists of one *Python* file with several classes providing the base functionality. As mentioned above, it relies on ADVACAM libraries for setting-up and reading the sensor. Other dependencies are well-known libraries from the “Python” eco-system for data analysis:

- `numpy`
- `matplotlib`,
- `scipy.ndimage.label`
- `numpy.cov`
- `numpy.linalg.eig`

The classes and scripts of the package are

- `class miniPIXdaq`
- `class frameAnalyzer`
- `class miniPIXana`
- `class runDAQ`
- `class bhist`
- `class scatterplot`
- `package script run_mPIXdaq.py`

Details on the interfaces are given below.

```
class miniPIXdaq:
    """Initialize, readout miniPIX device and store data

    After initialization, data from the device is stored in a
    ring buffer and the current buffer index is sent to the
    calling process via a Queue in an infinite loop, which
    ends when data is entered in a command Queue.

    Args:

        - ac_count: number of frames to overlay
        - ac_time: acquisition time

    Queues for communication and synchronization

        - dataQ: Queue to transfer data
        - cmsQ: command Queue

    Data structure:

        - fBuffer: ring buffer with recent frame data

    """
class frameAnalyzer:
    """Analyze frame data
    - find clusters
    - compute cluster energies
    - compute position and covariance matrix of x- and y-coordinates
    - analyze cluster shape (using eigenvalues of covariance matrix)
    - construct a tuple with cluster properties
    - optionally write cluster data to a file in csv format
```

Note: this algorithm only works if clusters do not overlap!

Args of `__call__()` method: a 2d-frame from the miniPIX

Returns:

- `n_pixels`: number of pixels with energy > 0
- `n_clusters`: number of clusters
- `n_cpixels`: number of pixels per cluster
- `circularity`: circularity per cluster (0. for linear, 1. for circular)
- `cluster_energies`: energy per cluster
- `single_energies`: energies in single pixels
- `self.clusters` is a tuple with properties per cluster with mean of x and y coordinates, number of pixels, energy, eigenvalues of covariance matrix and orientation ($[-\pi/2, \pi/2]$) and the minimal and maximal eigenvalues of the covariance matrix of the energy distribution:

format of the tuple:

```
        ( (x, y), n_pix, energy, (var_mx, var_mn), angle, (xE_m, yE_m), (varE_mx, varE_mn) )
"""
```

```
class miniPIXana:
```

```
    """Analysis of miniPIX frames for low-rate scenarios,
    where on-line analysis is possible and animated graphs are meaningful
```

```
    Animated graph of (overlaid) pixel images and cluster properties
```

```
    """
```

```
    Args:
```

- `npix`: number of pixels per axis (256)
- `nover`: number of frames to overlay
- `unit`: unit of energy measurement ("keV" or " μ s ToT")
- `circ`: circularity of "round" clusters (0. - 1.)
- `acq_time`: accumulation time per read-out frame

```
    """
```

Objects of these classes are instantiated by the class `runDAQ`. This class also accepts the command-line arguments to set various options, as already described above.

```
class runDAQ:
```

```
    """run miniPIX data acquisition and analysis
```

```
    class to handle:
```

- command-line arguments
- event loop controlling data acquisition and data output to file
- instantiates classes and calls corresponding methods for
 - initialization of miniPIX device
 - real-time analysis of data frames
 - animated figures to show a live view of incoming data

```
    """
```

Two helper classes implement 1d and 2d histogramming functionality for efficient and fast animation using methods from `matplotlib.pyplot`.

```
class bhist:
```

```
    """one-dimensional histogram for animation, based on bar graph
    supports multiple classes as stacked histogram
```

```
    Args:
```

```

* data: tuple of arrays to be histogrammed
* binedges: array of bin edges
* xlabel: label for x-axis
* ylabel: label for y axis
* yscale: "lin" or "log" scale
* labels: labels for classes
* colors: colors corresponding to labels
"""

```

```

class scatterplot:

```

```

    """two-dimensional scatter plot for animation, based on numpy.histogram2d
    supports multiple classes of data, plots a '.' in the corresponding color
    in every non-zero bin of a 2d-histogram

```

```

    Args:

```

```

        * data: tuple of pairs of coordinates (([x], [y]), ([], []), ...)
          per class to be shown
        * binedges: 2 arrays of bin edges ([bex], [bey])
        * xlabel: label for x-axis
        * ylabel: label for y axis
        * labels: labels for classes
        * colors: colors corresponding to labels
    """

```

A package script `run_mPIXdaq` is provided as an example to tie everything together in a running program. Because the ADVACAM *Python* interface (`pypixet.so`) expects C-libraries and configuration files in the very same directory as `pypixet.so` itself, some tricky manipulation of the environment variable `LD_LIBRARY_PATH` is needed to ensure that all libraries are loaded and the *miniPIX* is correctly initialized.

```

#!/usr/bin/env python3
import os, sys

# pypixet requires '.' in LD_LIBRARY_PATH so that all necessary C-libraries are found
# - add current directory to LD_LIBRARY_PATH
# - and restart python script for changes to take effect

path_modified = False

if 'LD_LIBRARY_PATH' not in os.environ:
    os.environ['LD_LIBRARY_PATH'] = '.'
    path_modified = True
elif not '.' in os.environ['LD_LIBRARY_PATH']:
    os.environ['LD_LIBRARY_PATH'] += ':'
    path_modified = True

if path_modified:
    print(" ! added '.' to LD_LIBRARY_PATH")
    try:
        os.execv(sys.argv[0], sys.argv)
    except Exception as e:
        sys.exit('EXCEPTION: Failed to Execute under modified environment, ' + e)
else:
    # restart python script for setting to take effect
    # get current working directory before importing minipix libraries
    wd = os.getcwd()
    from mpixdaq import mpixdaq # this changes the working directory!

    rD = mpixdaq.runDAQ(wd) # start daq in working directory
    rD()

```

Note that *Python* 3.7.9 is required to run under Microsoft Windows. Changing the environment does not work either under Windows, and therefore a simplified version of the run-script must be used:

```
from mpixdaq import mpixdaq  
rD = mpixdaq.runDAQ()  
rD()
```

It is also possible to start the script as a *Python* module:

```
python -m mpixdaq
```